


# ELECTROPHORESIS AUTOMATED SYSTEM

## ADALYA 40

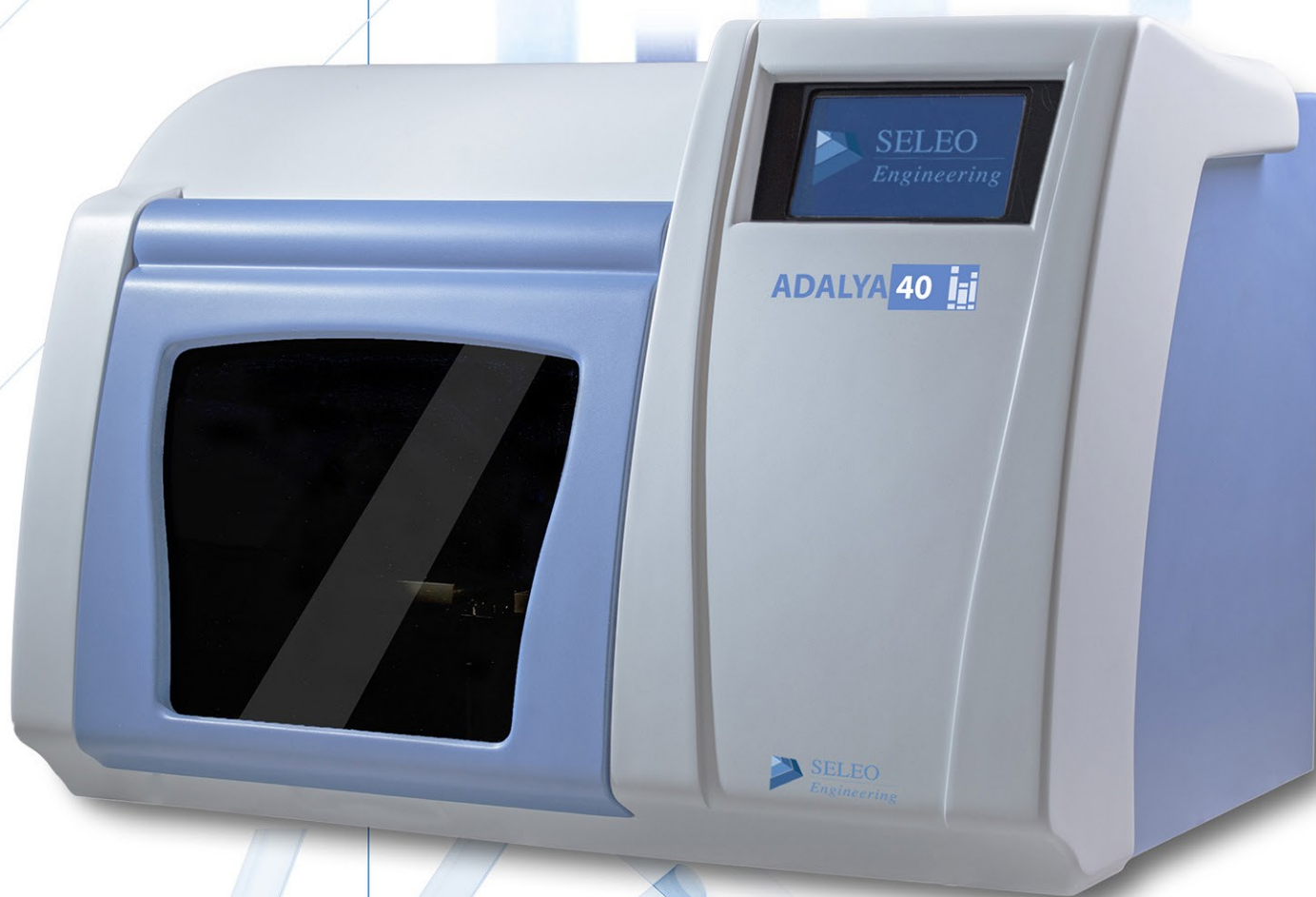


**SELEO**  
*Engineering*

Seleo Engineering s.r.l.  
1<sup>a</sup> traversa bugnano s.c.  
81030 Orta di Atella (CE) Italy  
TEL +39 (0)81 5022877 FAX +39(0)81 502287702  
Homepage [www.seleo.it](http://www.seleo.it)  
e-mail [info@seleo.it](mailto:info@seleo.it)

**MADE IN  
ITALY**

# ELECTROPHORESIS AUTOMATED SYSTEM

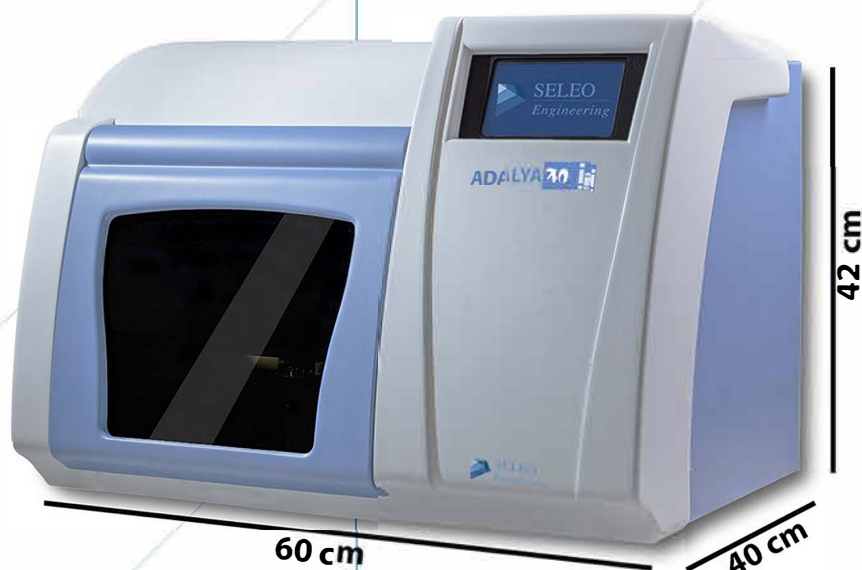

## FAST/VERSATILE

4/40 samples, up to 5  
electrophoresis membranes  
8 samples each.

Allows to perform the most common  
techniques:

- Serum Proteins
- Haemoglobins
- Lipoproteins
- Urinary Proteins

## COMPLETE AUTOMATION

Adalya automatically arranges  
the loading of all the reagents  
to standardize the delicate  
electrophoretic process.

During test all the liquid are  
constantly measured and  
adjusted, minimising the  
variables introduced by reagents  
level variations and  
to ensure the test preparation in true  
complete automation.



Serum Proteins



Haemoglobins



Lipoproteins

## PRODUCTIVITY

Serum-Proteins : 24 samples/ 1h:15'

Lipoproteins: 16 samples/ 1h:15'

Haemoglobins: 16 samples/ 1h:15'

## SOFTWARE FEATURES

User interface: monitor, keyboard and mouse

Minima recognition: automatic

Managing techniques: adjustable techniques in times, migration voltage/current

Operating system:

